

# The ultimate mobility tool

**Batec handbikes** are a comprehensive solution to your mobility problems:

- ✓ Go from reduced to **augmented mobility**; from needing help to being **independent** and **proactive**.
- ✓ Get a **2 in 1**, without giving up your manual wheelchair.
- ✓ **Overcome obstacles** and **drive over bumpy terrains**.
- ✓ Attach our handbikes **in just seconds** and **with no need for transfers**.
- ✓ Also available for **quadriplegics** and **hemiplegics**.



Request a quote or a demonstration  
without any commitment!



You CAN get used  
to being in a  
**wheelchair**



But not to giving up on

**LIFE**

[www.batec-mobility.com](http://www.batec-mobility.com)

**BATEC MOBILITY S.L.**

Illa de Buda, 2

08192 Sant Quirze del Vallès

(Barcelona) - ESPAÑA

93 544 2003

[info@batec-mobility.com](mailto:info@batec-mobility.com)

**BATEC**  
MOBILITY



# BATEC SCRAMBLER 2

Adventure is everywhere



# BATEC SCRAMBLER 2

## Adventure is everywhere

More powerful and safer, the **BATEC SCRAMBLER 2** arrives with improvements that will delight **fans of speed and off-road experiences**. If you are looking for intense sensations, the SCRAMBLER 2 is the perfect option for you. With innovative proprietary electronics and an even more off-road look.




All-in-one LCD screen with gear selector. Horn and USB port.



New proprietary electronics and anti-slip function.



300 RPM 1440 W high-torque brushless motor with reverse gear. Maximum speed 30 km/h.

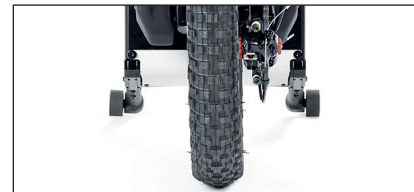

19" wheel with ultra-grip tyre and double-wall black aluminium rim.



200/203 mm disc, two mechanical brakes, two aluminium levers with parking brake.



Stand frame with suspension. Includes double LED taillight.



2000 lumens dual frontlight. Position and high beam lights. Retro-inspired yellow lens.



Storage bag.



55 km\* autonomy with 48 V – 768 Wh battery.

Battery with padded protective stickers in a matte black finish



Matte black aluminium handlebar with folding system for convenience.




G&B handlebar option without levers for quadriplegics.



Handlebar option with right or left Monolever brakes for hemiplegics.

### MEASUREMENTS

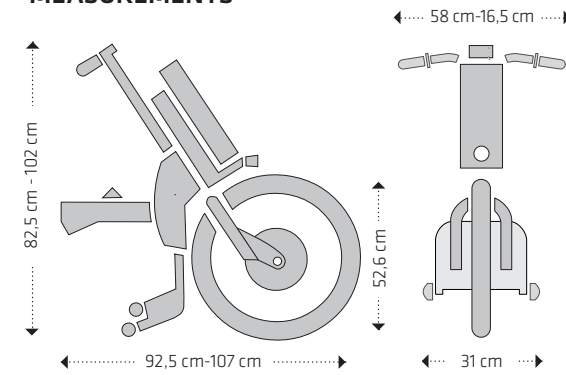

### COLOUR



\* Autonomy calculated under the following conditions: 75 - kg user, ultra-light wheelchair, flat terrain, no wind and steady speed of 15 km/h. May vary +/- 20% depending on the battery.  
Number/version: MM-009/4.0 | February 2023